

DETRO[®] AIR SPRAY

Innovative Technology

Airborne Disinfection

High level
disinfection
minimum chemical
usage



www.detrox.com.tr



Detrox Turkey
Atatürk mah. Adnan Menderes cad. No:7 Esenyurt /Istanbul/Turkey
Tel: +90 (212) 659 77 62 Fax: +90 (212) 659 77 63
e-mail: support@detrox.com.tr - international@detrox.com.tr

Detrox GmbH
Innstr. 16 D-68199 Mannheim/Germany
Tel: +49 (0) 621 899 585 04 Fax: +49 (0) 621 899 585 03
e-mail: info@detroxgmbh.de web: www.detroxgmbh.de



www.detrox.com.tr

High level disinfection
minimum chemical usage




DETRO[®]AIR SPRAY

PRODUCT FEATURES

- It's proper for disinfection of areas between 1 - 1.000 m³.
- Easy to start and adjust the program to be applied through the keypad on.
- Auto printing action through it's integrated printer. It can also print 100 pcs of operation retrospectively.
- It's easy to operate the machine through the short cut keys and it gives the chance to the user to save records for up to 100 rooms through "List Recording" function.
- It gives the chance to set the timing for start up.
- It can be controlled remotely from to distance of 50-100 meters thorough the optional RF Controller.
- The operation can be paused and resumed just in case.
- Adjustable screen brightness.
- Passwords can be adjusted in order to authorize the staff by whom the machine can be used.
- In order to perform easy and effective pulverization the head is positioned by the angel of 45°.
- Visible and audible warning system at the begining and the end of the operation.
- Turkish & English Languages are loaded up as default.
- Easy to hold through the plastic grip.
- Particle filter keeps the air access system clean.
- ABS Body
- Metal Parts are made of 304 Quality Stainless Steel
- 3X1.5 mm- 2 mt Length Power Supply Cable

TECHNICAL FEATURES

- Voltage: 220-240 VAC 50/60 Hz
- Weight: 7 kg
- Power: 1100W ENGINE
- Dimensions: 360 X 466 X 450 (width X length X height)

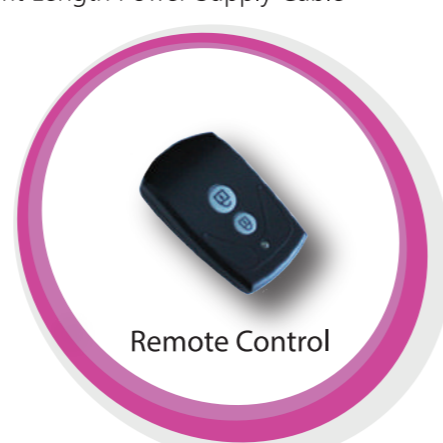

DETROSAN[®]AIR

PRODUCT SPECIFICATIONS

- Stabilized hydrogen peroxide based disinfectant solution for automated airborne room disinfection bactericidal and yeasticidal spectrum.
- Aldehyde free
- Leves no trace, wetness or gressy film
- %99 biodegradable, ecofriendly formula
- Non allergetic
- Bactericidal, fungicidal, virucidal, sporicidal
- Wide compatibility with surfaces and materials, noncorrosive

COMPOSITION

Stabilized hydrogen peroxide, silver salt complexes

INDICATIONS

Airborne disinfection of pre cleaned surfaces (Terminal disinfection)

ACTION MECHANISM

- Hydrogen peroxide (H₂O₂) is converted quickly into hydroxyl radicals (OH⁻) those react immeditely with the essential biomolecules (proteins, lipids, nucleic acids, etc.) inactivating them.
- Silver cations (Ag⁺), binding particularly with proteins, denaturing enzymes and interfering with DNA replication.
- The combined action of two active principles leads to cell inactivation to destruction of viruses, bacterias and spores.

INSTRUCTIONS FOR USE

- Ready for use solution to be used with Detro Air Spray
- Should be used in enclosed rooms without human presence
- Do not rinse

PACKING

- 1 L, 5 L, 20 L



EFFICIENCY

	Standarts
Bactericidal	EN 13727
Fungicidal	EN 13624
Virucidal	EN 14476
Tuberculocidal	EN 14348
Sporicidal	EN 13704

